




HYDRO-BIOS

Smart Sampling

Datasheet V-Fin Depressor 22kg (438 416)

V-Fin Depressor 22kg

Inverse vane for horizontal plankton nets. The V-Fin Depressor gives negative lift to the towing rope during horizontal or oblique plankton net operation. This also stabilises the net in the water and allows for good sampling results. It consists of a lead-weighted body and a stainless steel rope with shackle.

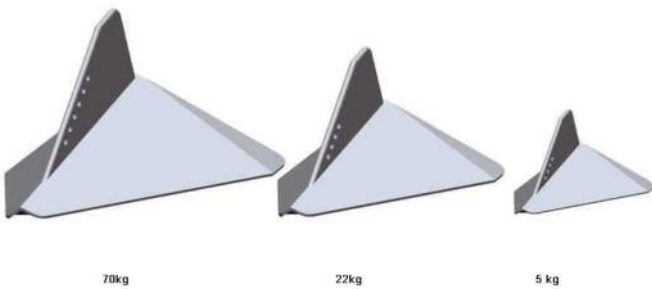

TECHNICAL DETAILS

length	620mm
width	670mm
height	400mm
weight on air	22kg
materials	anodized aluminium, lead
product ID	438 416